
BORT TaloStabil Plus Sport Aktiv-Sprunggelenkbandage

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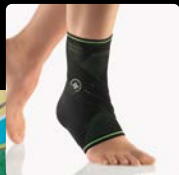Sport


Illustrationen
Figures

1


2

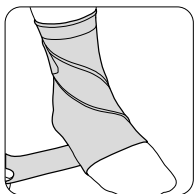
3

4

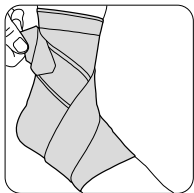
5

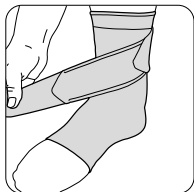6

7

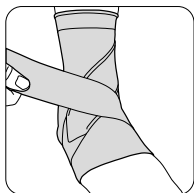8

Sprachen Languages

DE	Deutsch	Gebrauchsanweisung	04
EN	English	Instructions for use	09
FR	Français	Mode d'emploi	14
ES	Español	Instrucciones de uso	19
IT	Italiano	Istruzioni per l'uso	24
NL	Nederlands	Gebruiksaanwijzing	29
CS	Čeština	Návod k použití	34
ET	Eesti	Kasutusjuhend	38
PL	Polski	Instrukcja użytkowania	42
RO	Română	Instrucțiunile de utilizare	47

PDF: ga.bort.com

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Many thanks for placing your trust in a medical device from BORT GmbH. Please read the existing instructions for use carefully. If you have any questions, please refer to your physician or the specialist retailer from whom you purchased this medical device.

Intended purpose

This medical device is a support for soft ankle tissue compression consisting of elastic material with additional functional elements. They serve to treat the stabilisation of the upper and lower ankle joint.

Indications

Chronic post-traumatic or post-operative soft tissue irritations in the ankle area, arthrosis, ligament weakness, after removal of plaster cast, ligament injuries, distortion of upper ankle joint, supination prophyllaxis.

Contraindications

Thrombosis risk, extreme varicosis, peripheral arterial occlusive disease (PAOD), lymph drainage disorders, also unclear soft tissue swellings distal to the aid positioned, sensory and circulatory disorders with the region of the body treated, skin diseases in the part of the body treated.

Application risks/Important notes ⚠

This medical device is a prescribable product. Discuss the use and duration with your treating physician. The expert staff from whom you have received the medical device will select the appropriate size and instruct you regarding its use.

- remove the medical device prior to radiological examinations
- remove it during longer periods of rest, e.g. when sleeping
- in case of numbness, loosen the medical device and remove it if necessary
- in case of persistent complaints, consult the physician or specialist retailer
- use the medical device in accordance with therapeutic needs
- only use other products simultaneously after consultation with your physician
- do not make any changes to the product

- do not wear it on open wounds
- do not use in case of intolerance of one of the materials used
- no re-use. This medical aid is intended for treating one patient only
- whilst wearing the support, please neither use any creams nor ointments on or around it as they can destroy the material.
- In order to avoid damages to the knitted support fabric, only attach the Stabilo strap after successfully determining the size and trial fitting. Goods damaged by the Velcro fastening cannot be returned.

Fitting instructions

Grasp the upper edge of the support and pull it over the foot until the tops of the ankles are in the rounded pad recesses.

Two options are available for the attachment of the Stabilo strap provided for additional stabilisation:

For stabilisation with a raised outer edge of the foot, e. g. after twisting the foot outwards. Viewed from the rear, the strap points outwards.

- 1** Guide the Stabilo strap on the outer side of the ankle joint forwards diagonally over the metatarsal downwards to the inner side of the foot.
- 2** In the next step, guide the strap through underneath the sole of the foot to the outer side of the foot.
- 3** Then, draw the Stabilo strap upwards from the outer side diagonally over the metatarsal to the inner side of the ankle joint again. Raise the outer edge of the foot slightly here.

- 4** Finally, fasten it onto the outer side from behind.
Tip: The tighter you tension the Stabilo strap, the better the stabilising effect.

For stabilisation with a raised inner edge of the foot, e. g. after twisting the foot inwards. Viewed from the rear, the strap points inwards.

- 5** Guide the Stabilo strap on the outer side of the ankle joint forwards diagonally over the metatarsal downwards to the outer side of the foot.
- 6** In the next step, guide the strap through underneath the sole of the foot to the inner side of the foot.
- 7** Then, draw the Stabilo strap upwards from the inside diagonally over the metatarsal to the outer side of the ankle joint again. Raise the inner edge of the foot slightly here.
- 8** Finally, fasten it onto the inner side from behind.

Tip: The tighter you tense the Stabilo strap, the greater the stabilising effect is.

Removal

Start by loosening the Stabilo strap guided around the foot. Then, grasp the lower edge of the support and take it off by pulling it downwards over the toes.

Tip: Close the Velcro fasteners in order to avoid early wear to the Velcro surface.

Material composition

Polyamide/COOLMAX® (PA), polyamide (PA), elastane/LYCRA® (EL), silicone (SI).
The sewn in textile label provides the precise material composition.

Cleaning information

 Delicate wash 30 °C  Do not bleach  Do not dry in a tumble dryer

 Do not iron  Do not clean chemically

Do not use fabric conditioner. Stretch back into shape and dry in the fresh air.
Close the Velcro fastenings to avoid damage to other laundry items.

Guarantee

The legal regulations of the country in which you acquired the product apply to the purchased product. Please contact your specialist retailer if you suspect a warranty claim. Please clean the product before submitting a warranty claim. If the enclosed instructions for use have not been properly observed, the warranty may be impaired or cancelled. The warranty does not cover use of the product inappropriate for the indication, non-observance of application risks, instructions and unauthorised modifications to the product.

Useful life/Lifetime of the product

The lifetime of the medical device is determined by its natural wear and tear if treated appropriately and as recommended.

Duty of notification

If a serious worsening of the state of health occurs when using the medical device, you are obliged to notify the specialist dealer or us as the manufacturer and the Medicines & Healthcare products Regulatory Agency (MHRA).

You can find our contact information in these instructions for use. You can find the contact information for the national authority responsible for your country under the following link: www.bort.com/md-eu-contact.

Disposal

Upon the termination of use, the product must be disposed of in accordance with the corresponding local requirements.

Declaration of conformity

We confirm that this device conforms with the requirements of REGULATION (EU) 2017/745 OF THE EUROPEAN PARLIAMENT AND OF THE COUNCIL. You can find the current declaration of conformity under the following link: www.bort.com/conformity

Status: 10.2023

MD Medical device |  Single patient – multiple use

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