# BORT MobiDig Flexion

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

**REF** 112 560

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D112560|2024-10|006 ML | Rev.01
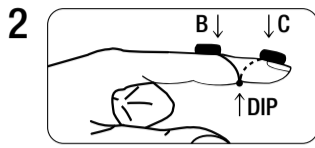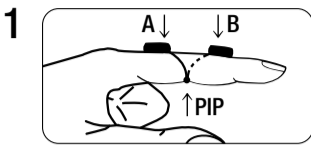
### DE BORT MobiDig Flexion

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Many thanks for placing your trust in a medical device from BORT GmbH. Please read the existing instructions for use carefully. If you have any questions, please refer to your physician or the specialist retailer from whom you purchased this medical device.

#### Intended purpose

This medical device is a finger brace for relieving the flexor tendon or for stimulating finger movement (dynamic Quengel casting).

#### Indications

Post-operative, post-traumatic flexion of the PIP or DIP joints.

#### Contraindications

Circulatory disorders with various origins, lymph drainage disorders, also unclear soft tissue swellings distal to the aid positioned, sensory disturbances of the region of the body treated, skin diseases in the part of the body treated.

#### Application risks/important notes

This medical device is a prescribable product. Discuss the use and duration with your treating physician. The expert staff from whom you have received the medical device will select the appropriate size and instruct you regarding its use.

- remove the medical device prior to radiological examinations
- remove it during longer periods of rest, e.g. when sleeping
- in case of numbness, loosen the medical device and remove it if necessary
- in case of persistent complaints, consult the physician or specialist retailer
- use the medical device in accordance with therapeutic needs
- only use other products simultaneously after consultation with your physician
- do not make any changes to the product
- do not wear it on open wounds
- do not use in case of intolerance of one of the materials used
- no re-use. This medical aid is intended for treating one patient only
- whilst wearing the support, please neither use any creams nor ointments on or around it as they can destroy the material

#### Fitting instructions

The finger splint/brace can be used for both the PIP joint (see Fig. 1) and the DIP joint (see Fig. 2). Tension can be increased or reduced by carefully bending the spring.

**1** Use as a PIP-Quengel: fit the splint between the proximal finger phalanx (A) and middle finger phalanx (B).

**2** Use as a DIP-Quengel: fit the splint between the middle finger phalanx (B) and end finger phalanx (C).

#### Removal

To remove the brace, turn it anticlockwise off the joint.

#### Material composition

Padding: 100% polyethylene foam (PE foam)  
Wire frame: spring steel type 1.4310

#### Cleaning information

Hand-washing 30 °C Do not bleach Do not dry in a tumble dryer  
 Do not iron Do not clean chemically  
Do not use fabric conditioner. If required, clean the brace with soap and lukewarm water. Air dry the brace.

#### Guarantee

The legal regulations of the country in which you acquired the product apply to the purchased product. Please contact your specialist retailer if you suspect a warranty claim. Please clean the product before submitting a warranty claim. If the enclosed instructions for use have not been properly observed, the warranty may be impaired or cancelled. The warranty does not cover use of the product inappropriate for the indication, non-observance of application risks, instructions and unauthorised modifications to the product.

#### Useful life/Lifetime of the product

The lifetime of the medical device is determined by its natural wear and tear if treated appropriately and as recommended.

#### Duty of notification

If a serious worsening of the state of health occurs when using the medical device, you are obliged to notify the specialist dealer or us as the manufacturer and the Medicines & Healthcare products Regulatory Agency (MHRA).

You can find our contact information in these instructions for use. You can find the contact information for the national authority responsible for your country under the following link: [www.bort.com/md-eu-contact](http://www.bort.com/md-eu-contact).

#### Disposal

Upon the termination of use, the product must be disposed of in accordance with the corresponding local requirements.

#### Declaration of conformity

We confirm that this device conforms with the requirements of REGULATION (EU) 2017/745 OF THE EUROPEAN PARLIAMENT AND OF THE COUNCIL. You can find the current declaration of conformity under the following link: [www.bort.com/conformity](http://www.bort.com/conformity)

Status: 10.2024

Medical device | Single patient – multiple use

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