

Novelties 2021



- ⊕ BORT StabiliGen **Eco**: The multi-variant solution in the series portfolio
- ⊕ Better fit, better function, more comfortable to wear

Med

In the midst
of healthy
life



With a new product in the
thumb category



NEW: Dynamic Osteoporosis
Brace



PediSoft: Now with new
Design

BORT. Benefit at your side.

Med

BORT StabiloGen Eco

The support for your standard treatment

INDIVIDUAL FIT

The even degree of extensibility of the knitted compression fabric adapts to the body shape.

FUNCTIONAL EFFECT

The compressive, massaging effect of the knitted fabric provides permanent proprioceptive, neuromuscular stimuli. This means that muscles are activated, the knee joint can be stimulated and pain relief be provided (gate control).

GOOD, PLEASANT WEARING COMFORT

Velvety surface due to the use of fluffy thread.

COMFORT ZONE POPLITEAL FOSSA

Highly elastic bending zone, which also protects the popliteal fossa when worn for longer periods.

VISCOELASTIC PAD

For stabilisation and relief of the knee joint

HIGH DIMENSIONAL STABILITY

One lateral spiral stay on each side

BORT. Benefit at your side.

Versatile functions and product variants

The multi-variant solution in the series portfolio

BORT StabiloGen Eco
REF 114 520



BORT StabiloGen Eco with a Grip Top Short
REF 114 520 HK



BORT StabiloGen Eco with a Grip Top
REF 114 520 H



Additional silicone grip top in all sizes,
height 5 cm



PRODUCT-BENEFIT

– For children too

- The classic product in standard format for the large majority of (straight) leg shapes
- Length of the support in size 3: 32 cm (12,6 inch)

COLOURS: skin-tone silver black denim

PRODUCT-BENEFIT

- Extensive thigh frame for short thigh lengths
- Additional silicone grip top, height 5 cm (1,97 inch), in all sizes
- Length of the support in size 3: 32 cm (12,6 inch)

COLOUR: silver

PRODUCT-BENEFIT

- Extensive thigh frame
- Additional silicone grip top, height 5 cm (1,97 inch), in all sizes
- Length of the support in size 3: 37 cm (14,57 inch)

COLOUR: silver

The original-Plus-size system

Better fit, better function,
more comfortable to wear

BORT PLUS-SIZE SYSTEM FOR THE THIGH AREA

- 5 Plus-sizes – As intermediate sizes in f-measure from size 2+ (40 cm, 15,75 inch) to size 6+ (68 cm, 26,77 inch) for a considerably better fit in the thigh area.
- Better function and more comfortable to wear.

SILICONE ADHESIVE BORDER FROM 5 SIZE



Keep your knee fit!

10 selected exercises for beginners and advanced users to improve mobility and muscular stability.



Off you go under:
uebungen.bort.com



BORT. Benefit at your side.

DYNAMIC OSTEOPOROSIS BRACE

Adapts to the course of therapy:
from the acute to the rehabilitation phase

BACK

BORT Generation Dynamic Osteoporosis Brace

REF 215 880

NEW



Generation

PRODUCT-BENEFIT

- Activating and spine-straightening back brace
- **Back splint and limitation restraint, stable and dynamically adjustable**
- Two-part, stable aluminium back splint, anatomically pre-shaped and individually adjustable without tools
- Relief and stabilisation of the lumbar spine
- Active straightening by activating the back muscles
- Increasing mobility

COLOUR: ■ blue

THE GENERATION-BENEFIT

- **Strap system according to the rucksack principle**
- Guide loops in the chest area for easy, stepless straightening of the upper body
- Lightweight and slim structure for good wearing comfort
- Fastener with a large finger hold and hand straps for effortless fitting

SHOULDER STRAP WITH REMOVABLE LIMITATION RESTRAINTS

- Fixed setting via the limitation restraints causes minimisation of the scope of movement (acute phase) and therefore supports a more stable straightening.
- The removal of the limitation restraints enables the realisation of dynamic reclination during the rehabilitation phase. Elastic shoulder strap parts effect more flexibility in the patient's scope of movement.

PROVEN RUCKSACK PRINCIPLE

- Shoulder strap with guidance fasteners for reclination eases putting it on and removing it without expending much force.
- The strength for reclination (straightening) can be adapted without expending much force.




TWO-PART ALUMINIUM BACK SPLINT FOR STABILISATION AND ACTIVATION

- The therapy with the dynamic brace must find the respective correct combination of stabilisation and activation.
- The back splint is joined to the additional reinforcement splint using clip fasteners for therapy during the acute phase.
- The brace is provided with more freedom of movement in the rehabilitation phase via the removal of the reinforcement plate.

DORSAL REDIRECTION OF THE SHOULDER STRAPS

- Dorsal, elastically held shoulder straps for better rotational freedom of movement.

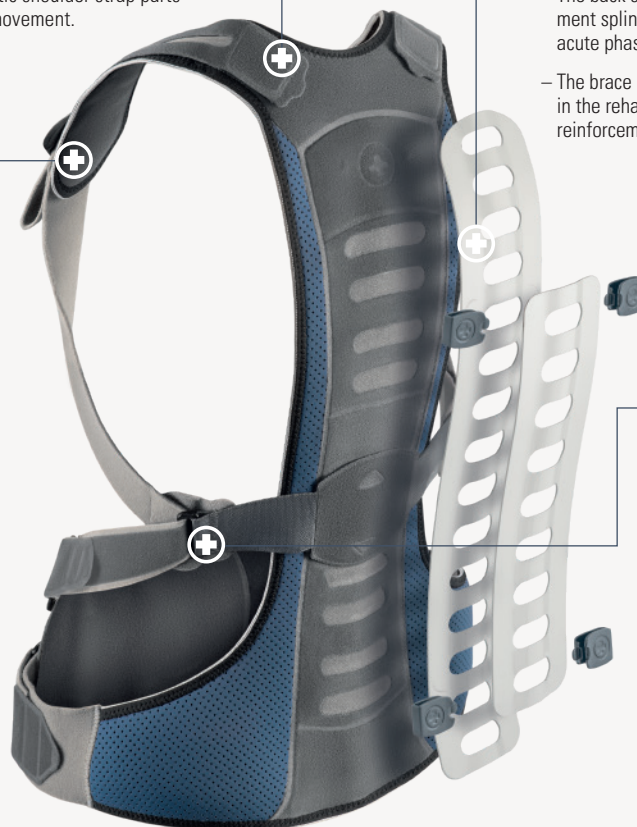

SOFT OA KNEE BRACE

For better mobility for the aged



KNEE



PRODUCT-BENEFIT

- Soft relief brace for isolated medial/lateral gonarthrosis axial difference up to 5° axial difference
- Adjustable, polycentric joint
- Limitation of extension: 0°, 5°, 10°, 15°, 20° and 25°
- Limitation of flexion: 0° (180°), 30°, 45°, 60°, 75° and 90°
- Brace length: 39 cm

COLOUR: ■ blue

THE GENERATION-BENEFIT

- Breathable material provides good wearing comfort
- Targeted relief of the medial/lateral compartment
- Varus/valgus correction freely adjustable
- Lightweight and slim structure for ideal wearing comfort

Correction of varus or valgus malalignment



Sport

THUMB SPLINT

The sporty side of BORT



THUMB



PRODUCT-BENEFIT

- Thumb brace for immobilisation of the metacarpophalangeal joint and saddle joint of the thumb
- Supportive in load situations

COLOUR: ■ black/green

THE SPORT-BENEFIT

- Breathable material provides optimal comfort
- Can be fully opened
- Can also be worn at night
- Modern, sporty design (with a reflecting label)



FINGER

BORT Buttonhole PIP Splint

REF 470 700

NEW



Med

PRODUCT-BENEFIT

- Plastic finger splint for treating buttonhole deformity (Boutonnière Deformity) on the long fingers of the hand
- Immobilisation of the PIP joint

COLOUR: ■ skin-tone

THE MED-BENEFIT

- Conservative treatment
- Thermoplastically mouldable
- Machining possible

Also available with the following units:

BORT Buttonhole PIP Splint 10 pack

REF 470 800

NEW

BORT Buttonhole PIP Splint set

REF 470 900

NEW

KNEE

BORT ROM Knee Brace

REF 100 900

NEW



Med

PRODUCT-BENEFIT

- Knee brace for mobilisation
- Brace length can be set using lateral telescopic splints
- Limitation of flexion:
0°, 10°, 20°, 30°, 40°, 50°, 60°, 70°, 80°, 90°, 100°, 110° and 120°
- Limitation of extension:
-10°, 0°, 10°, 20°, 30°, 40°, 50°, 60°, 70°, 80° and 90°
- Immobilisation:
-10°, 0°, 10°, 20°, 30°

COLOUR: ■ grey

THE MED-BENEFIT

- Setting joint limit stops and brace length without tools
- Lightweight and slim structure for ideal wearing comfort
- easyClick fasteners make it easy to put on/take off
- Easy to handle

PRODUCT MODIFICATIONS

SIZES




- ⊕ New size xx-large



- ⊕ New size xx-large




- ⊕ New size xx-large

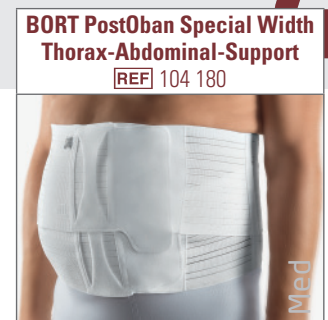

- ⊕ New size 0, with heights of 21 & 26 cm

CHANGED STRUCTURE/NEW DESIGN

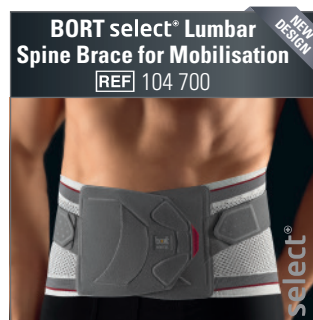

- ⊕ Easier to put on using the ergonomic latch plate



- ⊕ Easier to put on using the ergonomic latch plate



- ⊕ Easier to put on using the ergonomic latch plate



- ⊕ Easier to put on using the ergonomic latch plate

NEW DESIGN



- ⊕ Seamless comfort zone in the antecubital fossa




- ⊕ Breathable functional comfort zone in the popliteal fossa




- ⊕ Extended version, modified frame construction


COLOUR

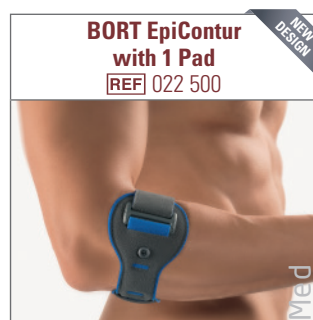

- ⊕ New colour grey
- tapered clasp design enables a non-slip, secure fit

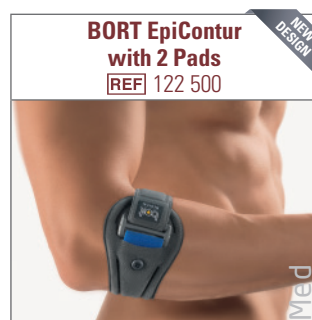

- ⊕ New colour grey
- tapered clasp design enables a non-slip, secure fit

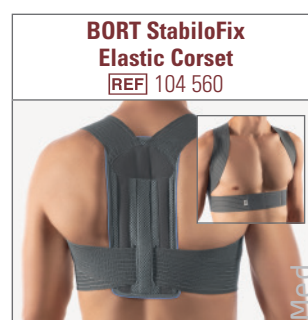

- ⊕ New colour grey

PEDISOFT

PEDI SOFT
Significantly better

NOW
WITH NEW
PACKAGING
DESIGN



We have revised our PediSoft lines for you and your customers under the motto "Significantly better"!

The modern packaging design can now score with its sustainable features.

The new PediSoft packaging consists of 93% recycled paper and is therefore compostable.

Your customers can simply dispose of the cardboard packaging, which is equipped with a film-free, natural water-based grease barrier, as waste paper.

Powerful POS presence:

The product range is rounded off with attractive displays in various sizes.



All information without guarantee, error and changes reserved. Status: January 2021.

BORT. Benefit at your side.

BORT GmbH P.O. Box 1330 | D-71367 Weinstadt | Am Schweizerbach 1 | D-71384 Weinstadt
Phone +49 (0) 7151 99200-0 | Fax +49 (0) 7151 99200-50 | www.bort.com | medical@bort.com

